

ProLong® Gold and *SlowFade*® Gold Antifade Reagents

Table 1. Contents and Storage Information.

Material	Amount	Storage	Stability
ProLong® Gold antifade reagent	1 bottle containing 10 mL (P36930)	<ul style="list-style-type: none"> • ≤-20°C • Protect from light 	When stored as directed, the product will be stable for at least 6 months
	5 bottles each containing 2 mL (P36934)		
ProLong® Gold antifade reagent with DAPI	1 bottle containing 10 mL (P36931)		
	5 bottles each containing 2 mL (P36935)		
<i>SlowFade</i> ® Gold antifade reagent	1 bottle containing 10 mL (S36936)		
	5 bottles each containing 2 mL (S36937)		
<i>SlowFade</i> ® Gold antifade reagent with DAPI	1 bottle containing 10 mL (S36938)		
	5 bottles each containing 2 mL (S36939)		

Introduction

Molecular Probes now offers improved versions of the very effective ProLong® and *SlowFade*® antifade reagents. ProLong® Gold antifade reagent (P36930, P36934) offers enhanced resistance to photobleaching and is premixed and ready to use—just add a drop of reagent and mount. ProLong® Gold reagent causes little or no quenching of the fluorescent signal (Figure 1). Like the original ProLong® antifade reagent, ProLong® Gold reagent will cure within 24 hours, and the sample can be saved for months after mounting. This antifade reagent offers excellent compatibility with a multitude of dyes and dye complexes (See *Appendix*, Table 2), making it an especially valuable tool for multicolor applications. For convenience, we also offer ProLong® Gold antifade reagent with DAPI (P36931, P36935); the addition of DAPI in the mounting medium eliminates the need for a separate nuclear counterstaining step.

SlowFade® Gold antifade reagent (S36936, S36937) outperforms the original *SlowFade*® reagent, offering increased resistance to photobleaching for a wide range of fluorescent dyes. *SlowFade*® Gold reagent is provided as a premixed and ready-to-use solution. Unlike ProLong® Gold reagent, *SlowFade*® Gold reagent does not cure over time, so samples can be viewed immediately—simply tack the corners of the slide with hot wax or nail polish, then image. *SlowFade*® Gold reagent is intended for short-term use (3–4 weeks) only; samples mounted using *SlowFade*® Gold reagent may degrade over time. *SlowFade*® Gold antifade reagent is also available with DAPI (S36938, S36939).

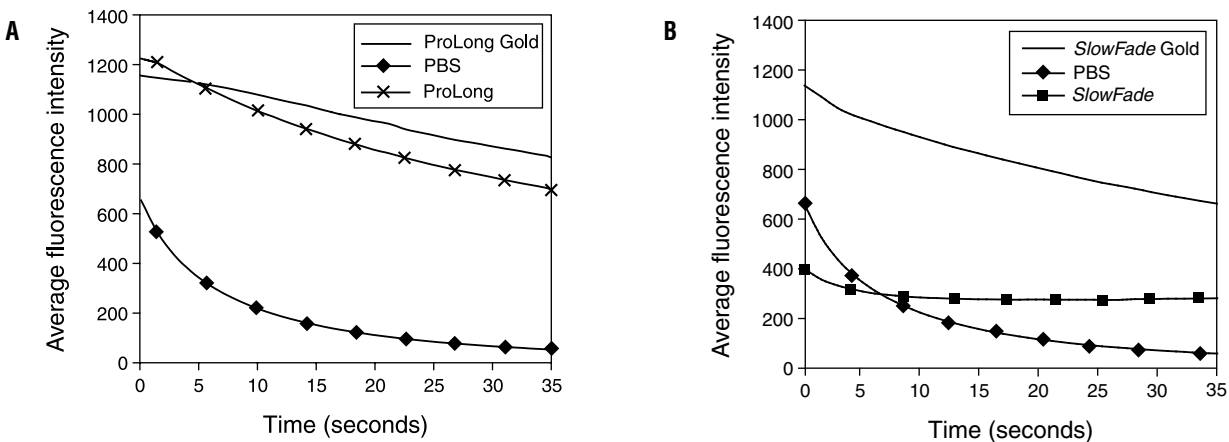


Figure 1. ProLong® Gold (Panel A) and SlowFade® Gold (Panel B) antifade reagents provide enhanced resistance to photobleaching. Fluorescein-labeled microspheres were mounted with various antifade reagents and illuminated for 30 seconds using a 100-watt Hg-arc lamp. Samples were imaged using a 40x/1.3 NA oil immersion lens and acquired using a 12-bit monochrome CCD camera. The same exposure settings were used for all images. The data plotted is the average fluorescence intensity from 20 microspheres.

Before You Begin

Important Considerations When Using ProLong® Gold Antifade Reagent

It is important to follow the instructions in this manual when using ProLong® Gold antifade reagent. The critical considerations for use of this product are summarized below:

- the bottle of ProLong® Gold antifade reagent **must** be warmed to room temperature before use
- excess moisture **must** be removed from the slide before ProLong® Gold antifade reagent is added
- after ProLong® Gold antifade reagent is added, the sample **must** be cured
 - place the mounted sample on a flat, dry surface
 - incubate for **24 hours at room temperature in the dark**

Viewing the Sample Briefly Before Curing

To view the sample briefly before curing, tack the corners of the coverslip with epoxy or nail polish. After viewing the sample, allow it to cure for 24 hours at room temperature in the dark on a flat, dry surface.

Extended Storage of Samples Mounted with ProLong® Gold Antifade Reagent

Following the curing time, the edges of the coverslip can be completely sealed with epoxy or nail polish, and the sample stored at room temperature, at 4°C, or at ≤-20°C. Sealing the edges retards the oxidation and extends the life of the sample for several months.

Technical Specifications for ProLong® Gold Antifade Reagent

- pH 7.4 at 20°C
- refractive index gradually increases as it cures (Figure 2)
- ProLong® Gold antifade reagent is useful for long-term storage (many months if edges are sealed), but must be cured for optimum performance

Important Considerations When Using *SlowFade*[®] Gold Antifade Reagent

It is important to follow the instructions in this manual when using *SlowFade*[®] Gold antifade reagent. The critical considerations for use of this product are summarized below:

- the bottle of *SlowFade*[®] Gold antifade reagent **must** be warmed to room temperature before use
- excess moisture **must** be removed from the slide before *SlowFade*[®] Gold antifade reagent is added

Viewing the Sample

Tack the corners of the coverslip with epoxy or nail polish, then image. *SlowFade*[®] Gold antifade reagent is intended for short-term use (3–4 weeks); samples mounted using this reagent may degrade over time

Technical Specifications for *SlowFade*[®] Gold Antifade Reagent

- pH 7.4 at 20°C
- refractive index is 1.42 at 20°C

Experimental Protocol

Protocol for Mounting Samples

These ready-to-use antifade reagents are ideal for use with most fixed samples. Because these solutions contain glycerol, they may be incompatible with some applications, such as mounting specimens that contain lipophilic plasma membrane stains like DiI.

1.1 Warm reagent. Remove the antifade reagent from the freezer and allow the vial to equilibrate to room temperature. Using an external heat source to warm the vial is not recommended, as this may decrease the long-term stability of the product.

1.2 Apply mounting medium. Remove any excess liquid from the specimen and apply 1 drop (or suitable quantity) of the antifade reagent to the specimen. Cover slide-mounted specimens with a coverslip; for specimens mounted on coverslips, place a drop of antifade reagent onto a clean slide and carefully lower the coverslip onto the antifade reagent to avoid trapping any air bubbles.

1.3 Prepare slides for viewing. For samples mounted using ProLong[®] Gold reagent, allow the preparation to cure on a flat surface in the dark. Curing time may vary from a couple of hours to overnight, depending on the thickness of the sample and the relative humidity of the surrounding air. We recommend 24 hours curing time. For long-term storage, seal the coverslip to the slide after curing to prevent excessive shrinkage of the mounting medium, which can result in sample distortion. After sealing, store the slide upright in a covered slide box at room temperature, at 4°C, or at ≤−20°C. Desiccant may be added to the box to ensure that the slide remains dry. To view the samples immediately, secure the coverslip at the corners using nail polish or hot wax to prevent the coverslip from moving. Leave the edges clear to allow the preparation to cure. After curing is complete, fully seal the edges for optimum sample longevity.

Samples mounted using *SlowFade*[®] Gold reagent can be imaged immediately after mounting. For thicker samples, it may take time for the antifade reagent to penetrate completely. If necessary, secure the coverslip at the corners using nail polish or hot wax to prevent the coverslip from moving.

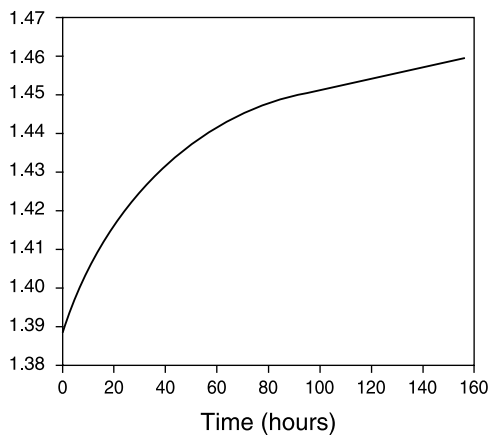


Figure 2. Increase in the refractive index of ProLong® Gold antifade reagent during the curing process.

Fluorescence Microscopy

Samples may be examined with a fluorescence microscope before the mounting medium dries. However, the antifade properties of ProLong® Gold and *SlowFade*® Gold antifade reagents do improve slightly the longer they are in contact with the specimen. ProLong® Gold reagent in particular achieves maximum effectiveness once it has cured. When properly stored, samples mounted in ProLong® Gold or *SlowFade*® Gold antifade reagent continue to resist photobleaching long after they are mounted.

To further reduce photobleaching, minimize the exposure of fluorescently labeled specimens to light by using neutral density filters, and expose samples only when observing or recording a signal. Optimize image capture by using a minimum of optics, high-numerical aperture objectives, relatively low magnification, high-quality optical filters, and high-speed film or high-efficiency detectors.

Appendix

Table 2. Resistance to photobleaching of dyes mounted in ProLong® Gold and SlowFade® Gold reagents.*

Dye	ProLong® Gold Reagent	SlowFade® Gold Reagent
Alexa Fluor® 350	65%	55%
Alexa Fluor® 488	85%	90%
Alexa Fluor® 546	67%	76%
Alexa Fluor® 555	80%	89%
Alexa Fluor® 568	93%	88%
Alexa Fluor® 594	92%	92%
Alexa Fluor® 633	99%	ND
Alexa Fluor® 635	100%	85%
Alexa Fluor® 647	100%	100%
Alexa Fluor® 660	94%	99%
Alexa Fluor® 680	96%	99%
7-Aminoactinomycin D	92%	91%
Cy2	89%	83%
Cy3	83%	83%
Cy5	100%	100%
Fluorescein	85%	78%
BODIPY® FL	54%	80%
Tetramethylrhodamine	98%	75%
Texas Red®	61%	75%
DAPI	67%	69%
SYTOX® Green	85%	77%
TO-PRO®-3	92%	72%

*Expressed as percentage of initial fluorescence intensity remaining following 30 seconds of continuous illumination using a 40X/1.3 NA objective and a 100-watt Hg-arc lamp as the light source. ND = Not determined.

Product List Current prices may be obtained from our website or from our Customer Service Department.

Cat #	Product Name	Unit Size
P36930	ProLong® Gold antifade reagent	10 mL
P36934	ProLong® Gold antifade reagent *special packaging*	5 X 2 mL
P36931	ProLong® Gold antifade reagent with DAPI	10 mL
P36935	ProLong® Gold antifade reagent with DAPI *special packaging*	5 X 2 mL
S36936	SlowFade® Gold antifade reagent	10 mL
S36937	SlowFade® Gold antifade reagent *special packaging*	5 X 2 mL
S36938	SlowFade® Gold antifade reagent with DAPI	10 mL
S36939	SlowFade® Gold antifade reagent with DAPI *special packaging*	5 X 2 mL

Contact Information

Molecular Probes, Inc.
29851 Willow Creek Road
Eugene, OR 97402
Phone: (541) 465-8300
Fax: (541) 335-0504

Customer Service:
6:00 am to 4:30 pm (Pacific Time)
Phone: (541) 335-0338
Fax: (541) 335-0305
probesorder@invitrogen.com

Toll-Free Ordering for USA:
Order Phone: (800) 438-2209
Order Fax: (800) 438-0228

Technical Service:
8:00 am to 4:00 pm (Pacific Time)
Phone: (541) 335-0353
Toll-Free (800) 438-2209
Fax: (541) 335-0238
probestech@invitrogen.com

Invitrogen European Headquarters
Invitrogen, Ltd.
3 Fountain Drive
Inchinnan Business Park
Paisley PA4 9RF, UK
Phone: +44 (0) 141 814 6100
Fax: +44 (0) 141 814 6260
Email: euroinfo@invitrogen.com
Technical Services: eurotech@invitrogen.com

Further information on Molecular Probes products, including product bibliographies, is available from your local distributor or directly from Molecular Probes. Customers in Europe, Africa and the Middle East should contact our office in Paisley, United Kingdom. All others should contact our Technical Assistance Department in Eugene, Oregon.

Molecular Probes products are high-quality reagents and materials intended for research purposes only. These products must be used by, or directly under the supervision of, a technically qualified individual experienced in handling potentially hazardous chemicals. Please read the Material Safety Data Sheet provided for each product; other regulatory considerations may apply.

Limited Use Label License

For research use only. Not intended for any animal or human therapeutic or diagnostic use. The purchase of this product conveys to the buyer the non-transferable right to use the purchased amount of the product and components of the product in research conducted by the buyer (whether the buyer is an academic or for-profit entity). The buyer cannot sell or otherwise transfer (a) this product (b) its components or (c) materials made using this product or its components to a third party or otherwise use this product or its components or materials made using this product or its components for Commercial Purposes. The buyer may transfer information or materials made through the use of this product to a scientific collaborator, provided that such transfer is not for any Commercial Purpose, and that such collaborator agrees in writing (a) to not transfer such materials to any third party, and (b) to use such transferred materials and/or information solely for research and not for Commercial Purposes. Commercial Purposes means any activity by a party for consideration and may include, but is not limited to: (1) use of the product or its components in manufacturing; (2) use of the product or its components to provide a service, information, or data; (3) use of the product or its components for therapeutic, diagnostic or prophylactic purposes; or (4) resale of the product or its components, whether or not such product or its components are resold for use in research. Invitrogen Corporation will not assert a claim against the buyer of infringement of the above patents based upon the manufacture, use or sale of a therapeutic, clinical diagnostic, vaccine or prophylactic product developed in research by the buyer in which this product or its components was employed, provided that neither this product nor any of its components was used in the manufacture of such product. If the purchaser is not willing to accept the limitations of this limited use statement, Invitrogen is willing to accept return of the product with a full refund. For information on purchasing a license to this product for purposes other than research, contact Molecular Probes, Inc., Business Development, 29851 Willow Creek Road, Eugene, OR 97402. Tel: (541) 465-8300. Fax: (541) 335-0504.

Several Molecular Probes products and product applications are covered by U.S. and foreign patents and patents pending. All names containing the designation ® are registered with the U.S. Patent and Trademark Office.

Copyright 2005, Molecular Probes, Inc. All rights reserved. This information is subject to change without notice.